



British Medical Auctions Ultrasound Equipment Quality Check (EQC) Preventive Maintenance Service Report (PM Service)

Esaote / MyLabX9 / 100012

Complete

British Medical Auctions

Address: Unit 1, Sterling Centre, Eastern Road, Bracknell, RG12 2PW
Tel: +44 1252 299024
Email: enquiries@britishmedicalauctions.co.uk
Web: www.britishmedicalauctions.co.uk

Department: Service

System has passed performance and safety assessments and is issued as fit for use, subject to final acceptance assessment upon fit.

Safe

System Type:

Ultrasound System



Photo 1

System Manufacture:

Esaote

Model:

MyLabX9

Serial No.:

100012



Photo 2

Year of Manufacture:

05.2021

Lot No.:

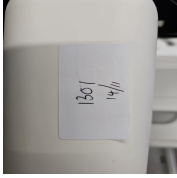


Photo 3

Software Level:

F120103

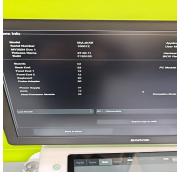


Photo 4



Photo 5

Service Code:

722

Electrical Safety Test Asset ID:

744

Test Equipment No. (Electrical Safety Analyzer):
Rigel 288+
SN.:10M-0263

Calibration Due Date:
05.2025

Test Equipment No. (Ultrasound Phantom):
ATS Laboratories MODEL: 539
Multipurpose Ultrasound Phantom
SN.: 200493148

Date Performed

19.11.2024 11:26 CET

Performed By

MS

1. PRELIMINARY CHECKS

Power Up

Yes



Photo 6

Database Cleared (Including all Patient Data and Previous Owner Details)

Yes

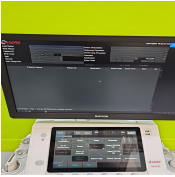


Photo 7

Password Protection

No

Time and Date Check

Yes

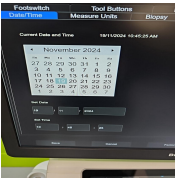


Photo 8

2. PHYSICAL INSPECTION

Input Power

Pass



Photo 9

External I/O

Pass

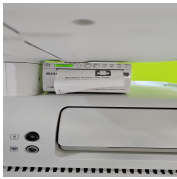


Photo 10

Wheels and Brakes

Pass



Photo 11



Photo 12



Photo 13



Photo 14

Control Panel Movement

Pass



Photo 15



Photo 16

Cables and Connectors

Pass

Shielding and Covers

Pass



Photo 17



Photo 18



Photo 19



Photo 20



Photo 21

Console

Pass



Photo 22



Photo 23

Control Panel

Pass



Photo 24

Monitor Arm

Pass



Photo 25

3. CLEANING

Monitor (CRT/LCD)

Yes



Photo 26

Touch screen

Yes

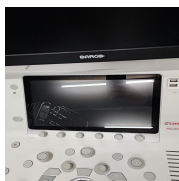


Photo 27

Console

Yes



Photo 28

Keyboard

Yes



Photo 29



Photo 30

Trackball / Optical mouse

Yes



Photo 31

Probe Holder

Yes



Photo 32



Photo 33

Internal Cleaning

Yes

Filter

Yes

Shielding and Covers

Yes



Photo 34



Photo 35



Photo 36

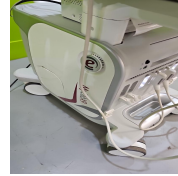


Photo 37



Photo 38

Thermal Printer Head

Yes

4. FUNCTIONAL CHECKS

Monitor

Pass



Photo 39



Photo 40

Touch Panel

Pass

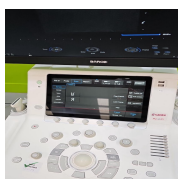


Photo 41

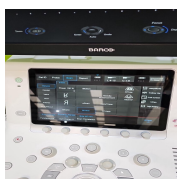


Photo 42

Console

Pass



Photo 43

Keyboard

Pass



Photo 44



Photo 45

B-Mode

Pass



Photo 46

M-Mode

Pass



Photo 47

C-Mode

Pass

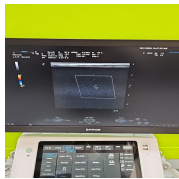


Photo 48

PD-Mode

Pass



Photo 49

PW-Mode

Pass

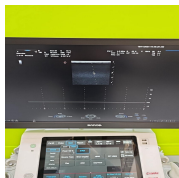


Photo 50

B/W Printer

Pass



Photo 51

5. PROBE NO.: 1

Probe Make

Esaote

Probe Model

C1-8

Probe Serial Number

2752



Photo 52

Visual inspection

Pass

Cable

Pass



Photo 53

Lens

Pass



Photo 54

Probe Casing

Pass



Photo 55



Photo 56



Photo 57



Photo 58

Probe Casing Strain Relief

Pass



Photo 59

Connector

Pass



Photo 60



Photo 61

Connector Strain Relief

Pass



Photo 62

Lateral resolution

Pass



Photo 63

Axial resolution

Pass



Photo 64

Image benchmark

Pass



Photo 65

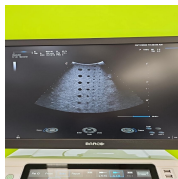


Photo 66

Probe passes inspection & test

Yes

Probe has passed performance and safety assessments and is issued as fit for use, subject to final acceptance assessment upon fit.

Yes

PROBE NO.: 2

Probe Make

Esaote

Probe Model

L4-15

Probe Serial Number

6685



Photo 67

Visual inspection

Pass

Cable

Pass



Photo 68

Lens

Pass



Photo 69

Probe Casing

Pass



Photo 70



Photo 71



Photo 72



Photo 73

Probe Casing Strain Relief

Pass



Photo 74

Connector

Pass



Photo 75



Photo 76

Connector Strain Relief

Pass



Photo 77

Lateral resolution

Pass



Photo 78

Axial resolution

Pass



Photo 79

Image benchmark

Pass



Photo 80



Photo 81

Probe passes inspection & test

Yes

Probe has passed performance and safety assessments and is issued as fit for use, subject to final acceptance assessment upon fit.

Yes

6. SUMMARY

This service has been completed in accordance with our quality procedure.
If you have any questions or concerns, please contact us on
+44 1252 299024
or email: enquiries@britishmedicalauctions.co.uk to discuss further.

System has passed performance and safety assessments and is issued as fit for use, subject to final acceptance assessment upon fit.

Safe

Date Performed:

19.11.2024 11:51 CET

Service Engineer Signature

MS
19.11.2024 11:51 CET

British Medical Auctions
Unit 1, Sterling Centre, Eastern Road, Bracknell, RG12 2PW
Tel: +44 1252 299024
Email: enquiries@britishmedicalauctions.co.uk
Web: www.britishmedicalauctions.co.uk

British Medical Auctions
Ultrasound Equipment Quality Check (EQC)
Preventive Maintenance Service Report (PM Service)

V1.0 05.2022

7. PHOTOGRAPHIC DOCUMENTATION



Photo 82



Photo 83



Photo 84



Photo 85



Photo 86

Media summary



Photo 1



Photo 2

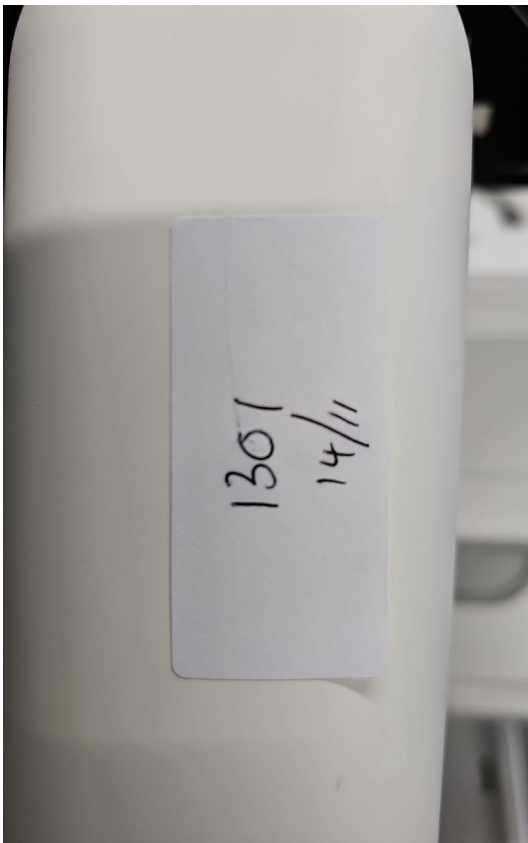


Photo 3

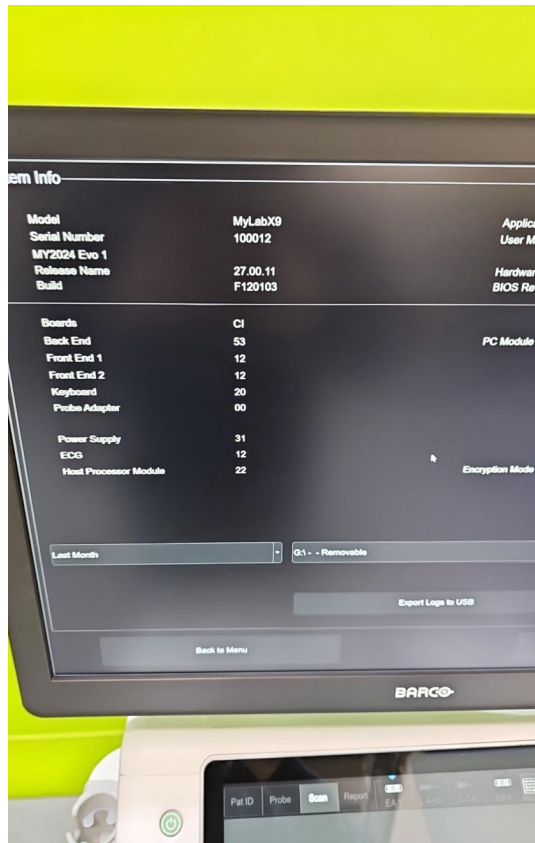


Photo 4

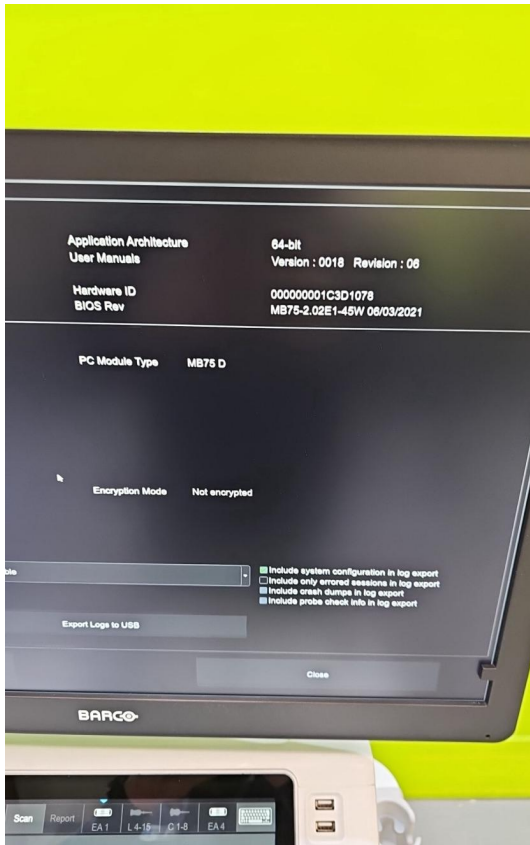


Photo 5



Photo 6

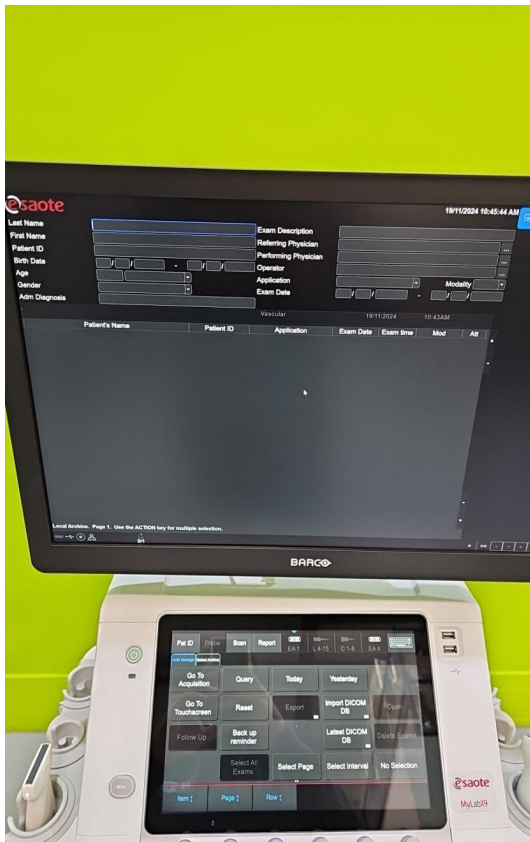


Photo 7

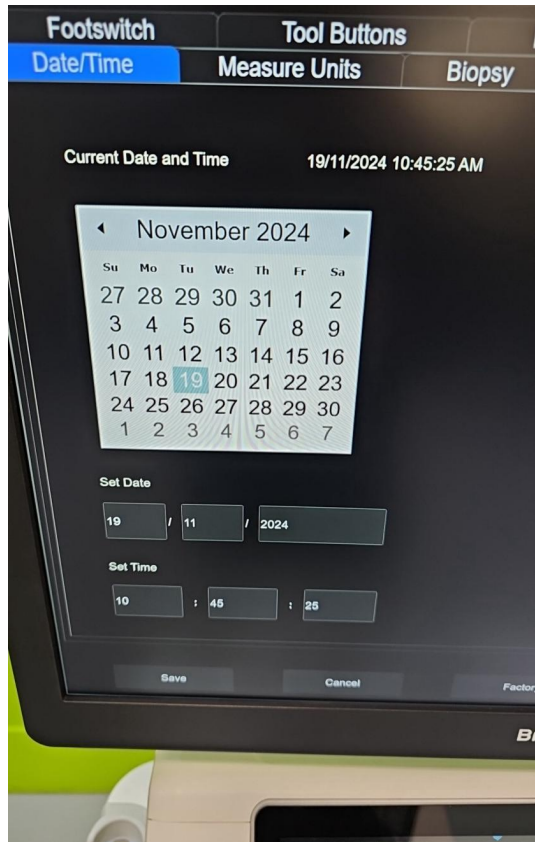


Photo 8



Photo 9

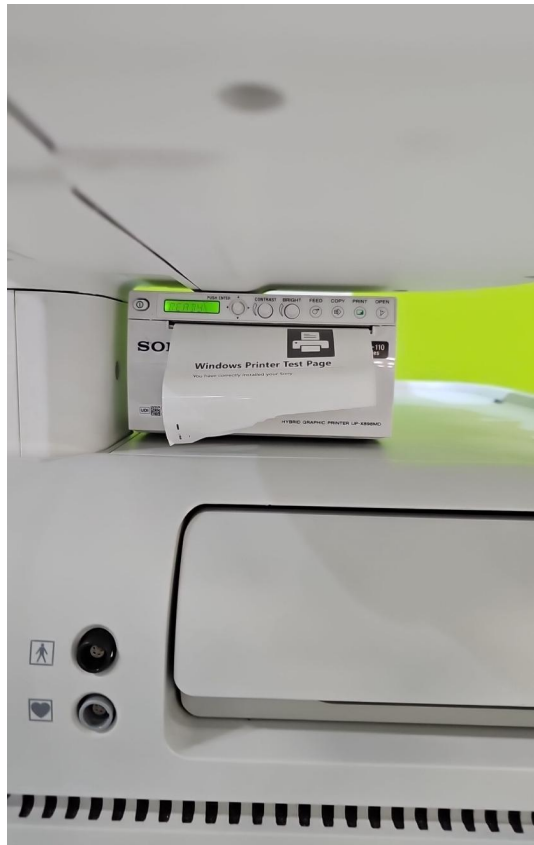


Photo 10



Photo 11



Photo 12



Photo 13



Photo 14



Photo 15



Photo 16



Photo 17



Photo 18



Photo 19



Photo 20



Photo 21



Photo 22



Photo 23



Photo 24



Photo 25

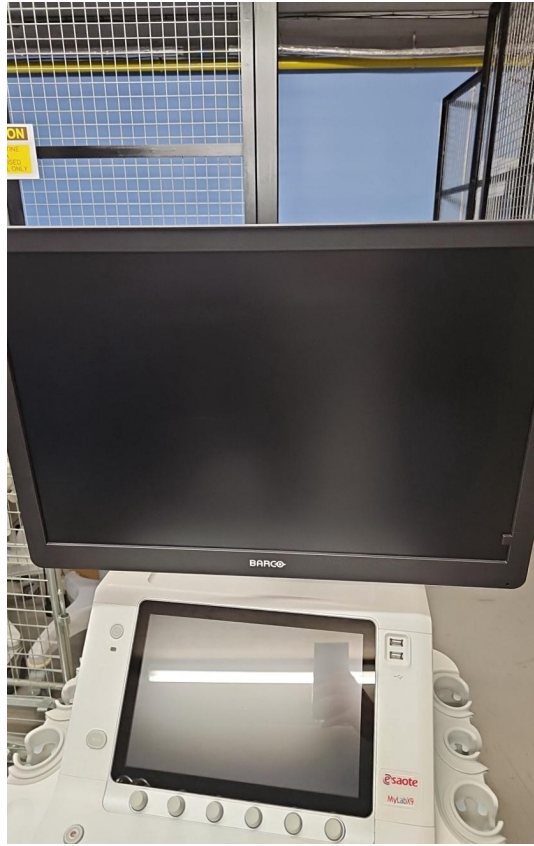


Photo 26

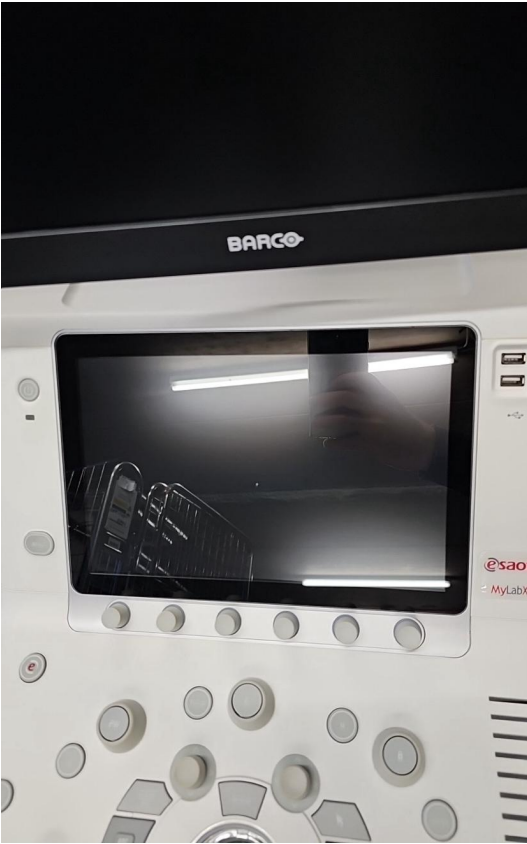


Photo 27



Photo 28

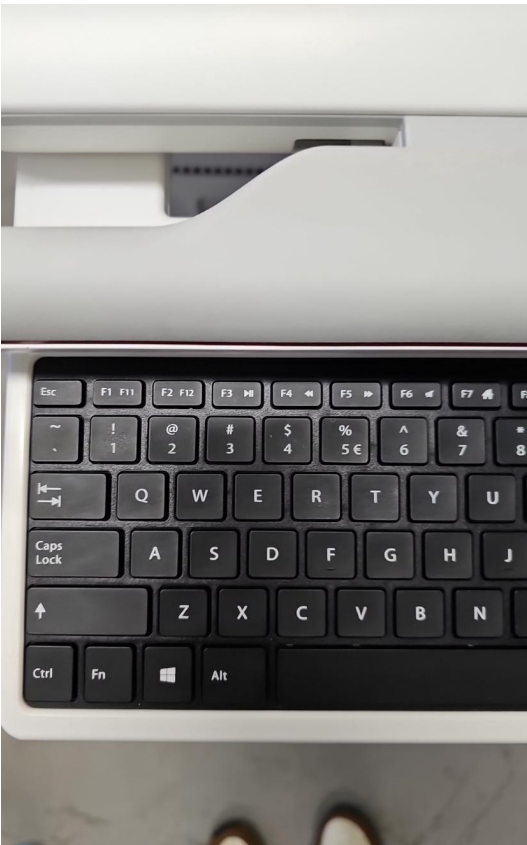


Photo 29

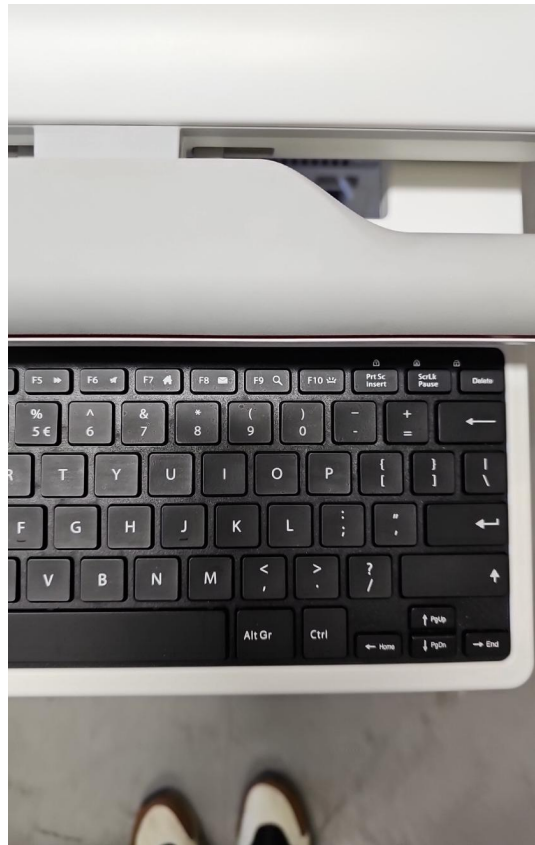


Photo 30



Photo 31



Photo 32



Photo 33



Photo 34



Photo 35



Photo 36



Photo 37



Photo 38



Photo 39



Photo 40



Photo 41

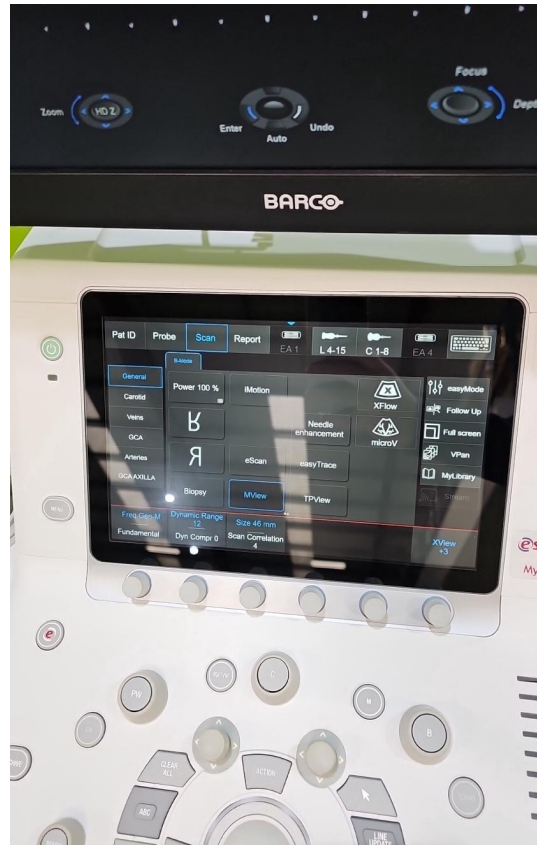


Photo 42

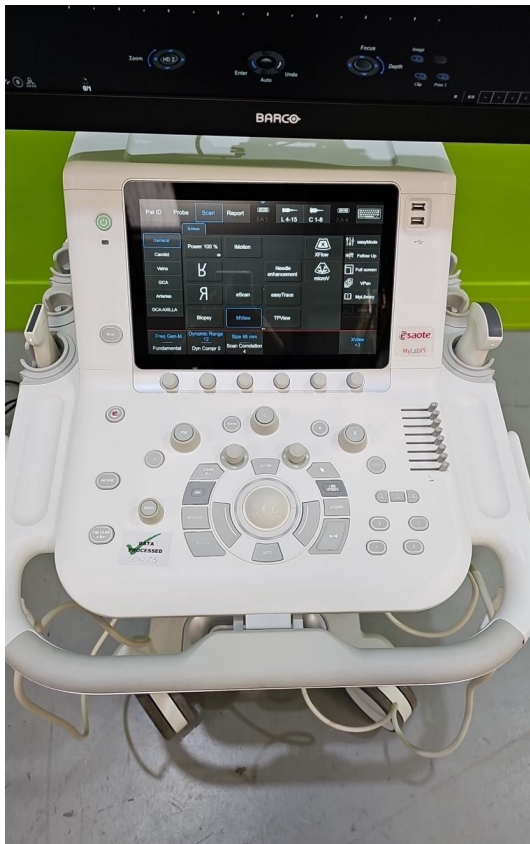


Photo 43



Photo 44



Photo 45

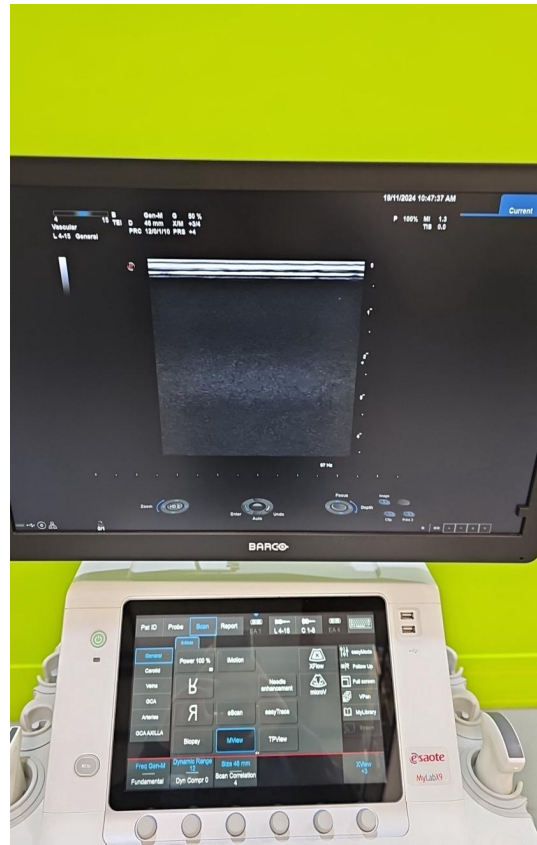


Photo 46

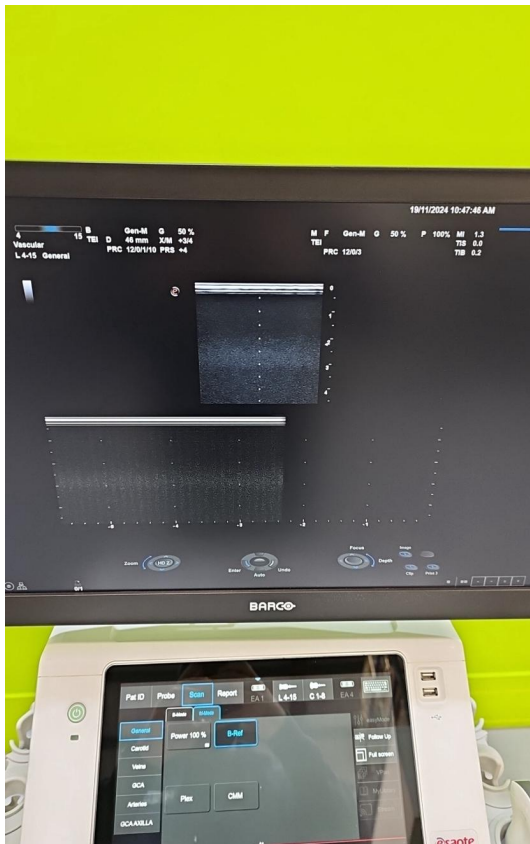


Photo 47

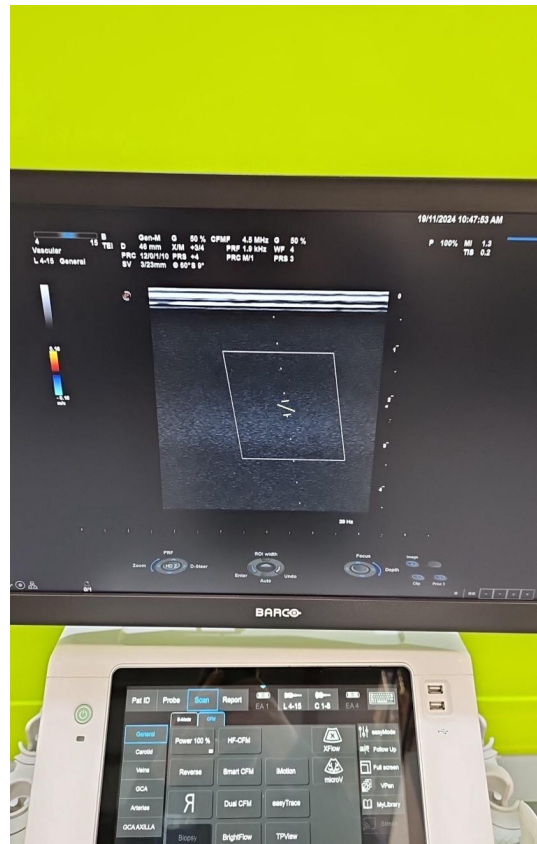


Photo 48

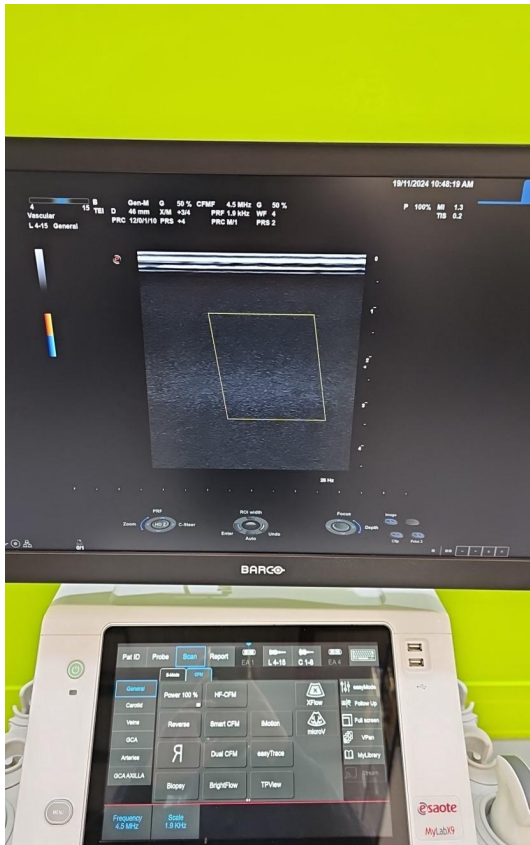


Photo 49

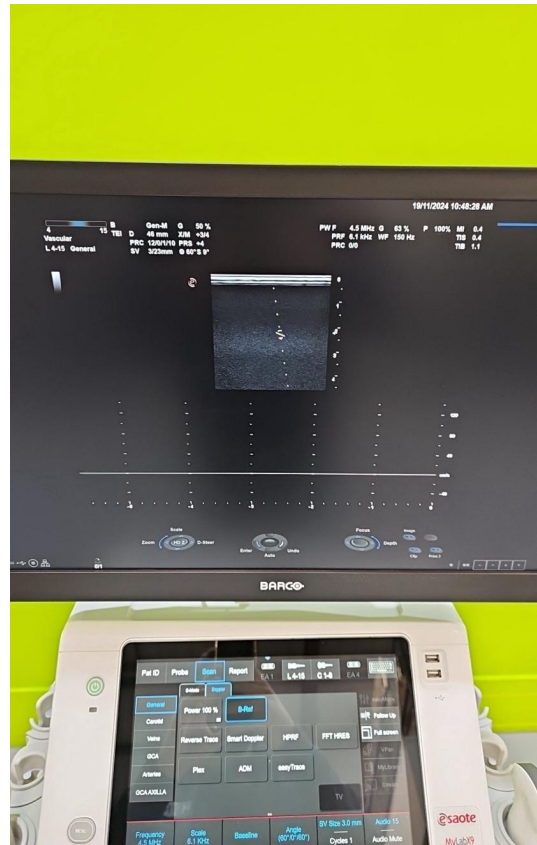


Photo 50

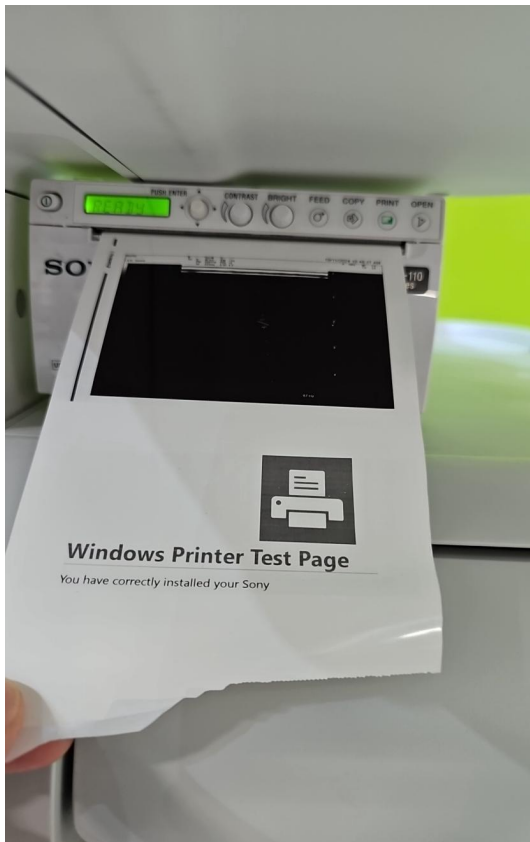


Photo 51



Photo 52



Photo 53



Photo 54



Photo 55



Photo 56



Photo 57

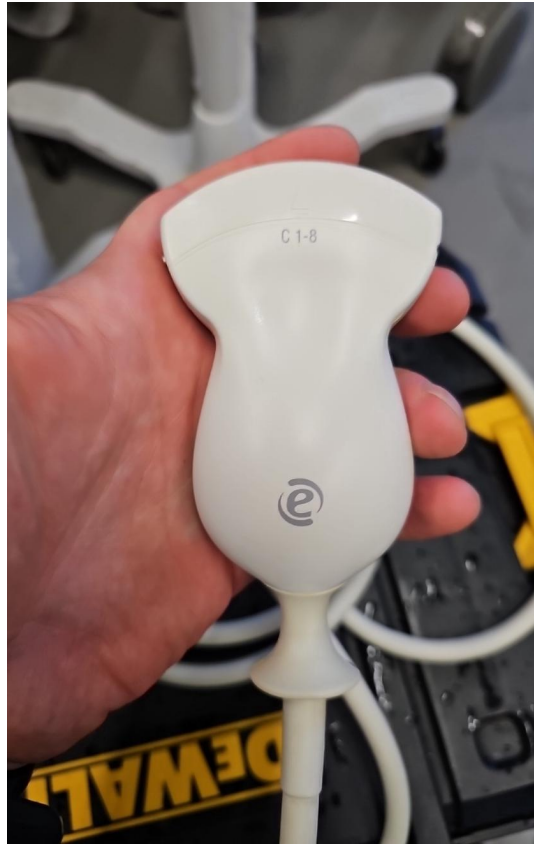


Photo 58



Photo 59



Photo 60



Photo 61



Photo 62



Photo 63



Photo 64

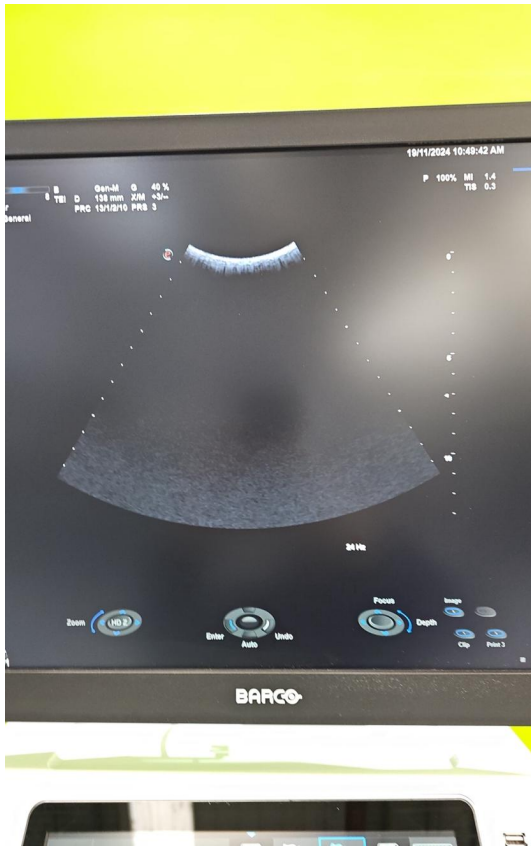


Photo 65

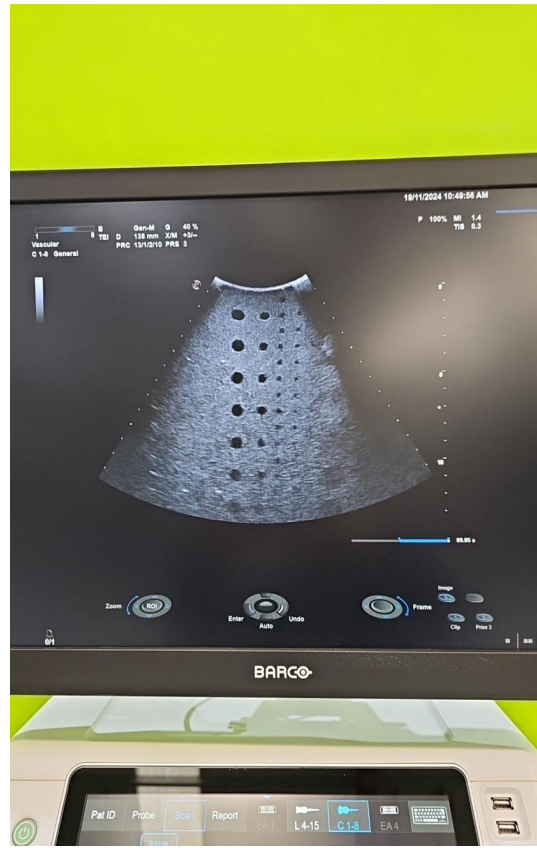


Photo 66



Photo 67



Photo 68



Photo 69



Photo 70



Photo 71



Photo 72



Photo 73



Photo 74



Photo 75



Photo 76



Photo 77

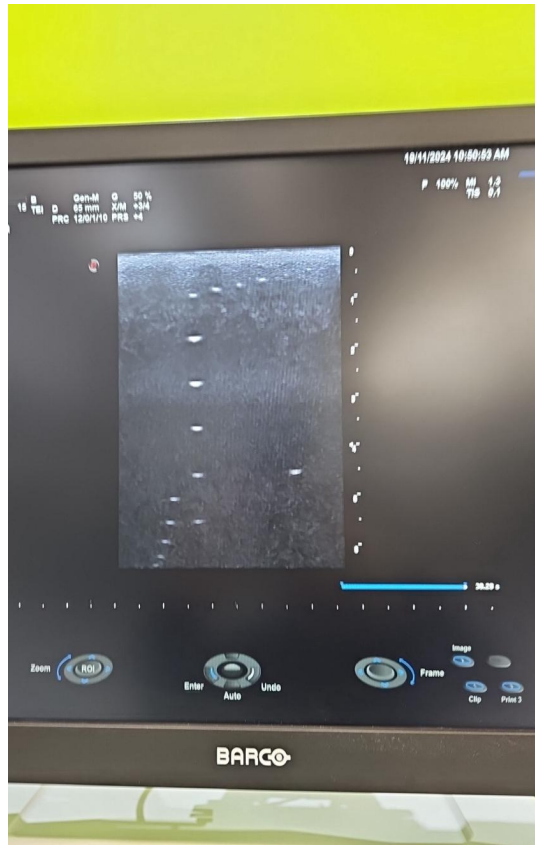


Photo 78

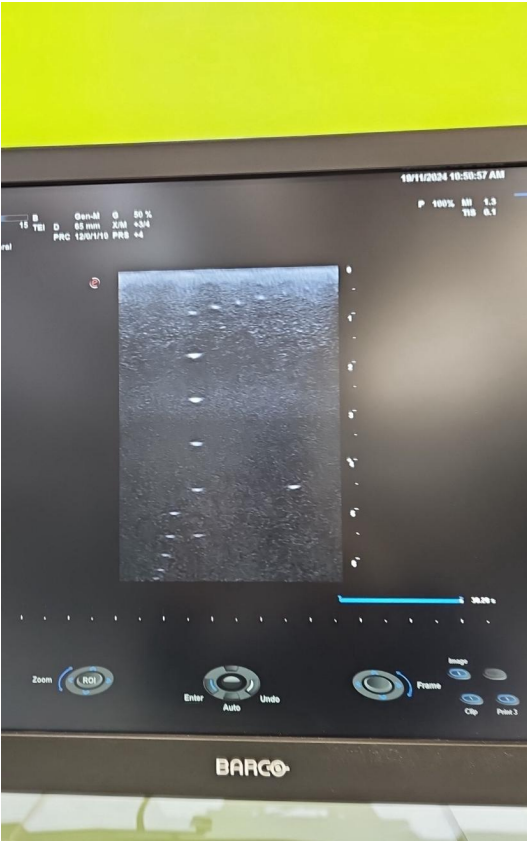


Photo 79

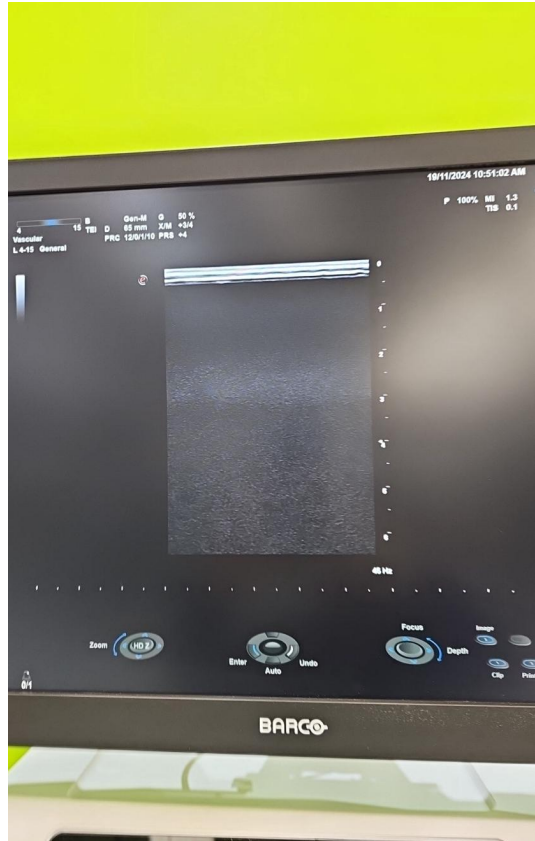


Photo 80



Photo 81



Photo 82



Photo 83



Photo 84



Photo 85

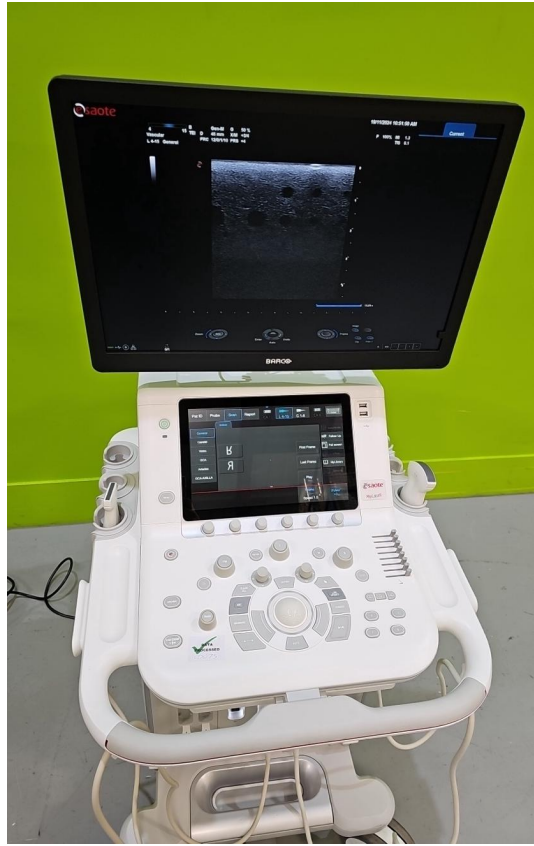


Photo 86